




Elite 3100A Top Slide Surgical Table




**Unobstructed C-Arm Access
for Full Body Imaging**

Maximum C-Arm Access with 23-1/2 inch Top Slide

Developed for maximum C-Arm access in procedures where only lateral tilt and Trendelenburg positioning are required, the Elite 3100A surgical table provides 23-1/2 inches of longitudinal top slide. This table is the ideal choice for vascular surgery, hip revisions, ureteroscopies, laser angioplasties, pacemaker insertions, pain management and other C-Arm assisted procedures that can be accomplished quickly and efficiently.

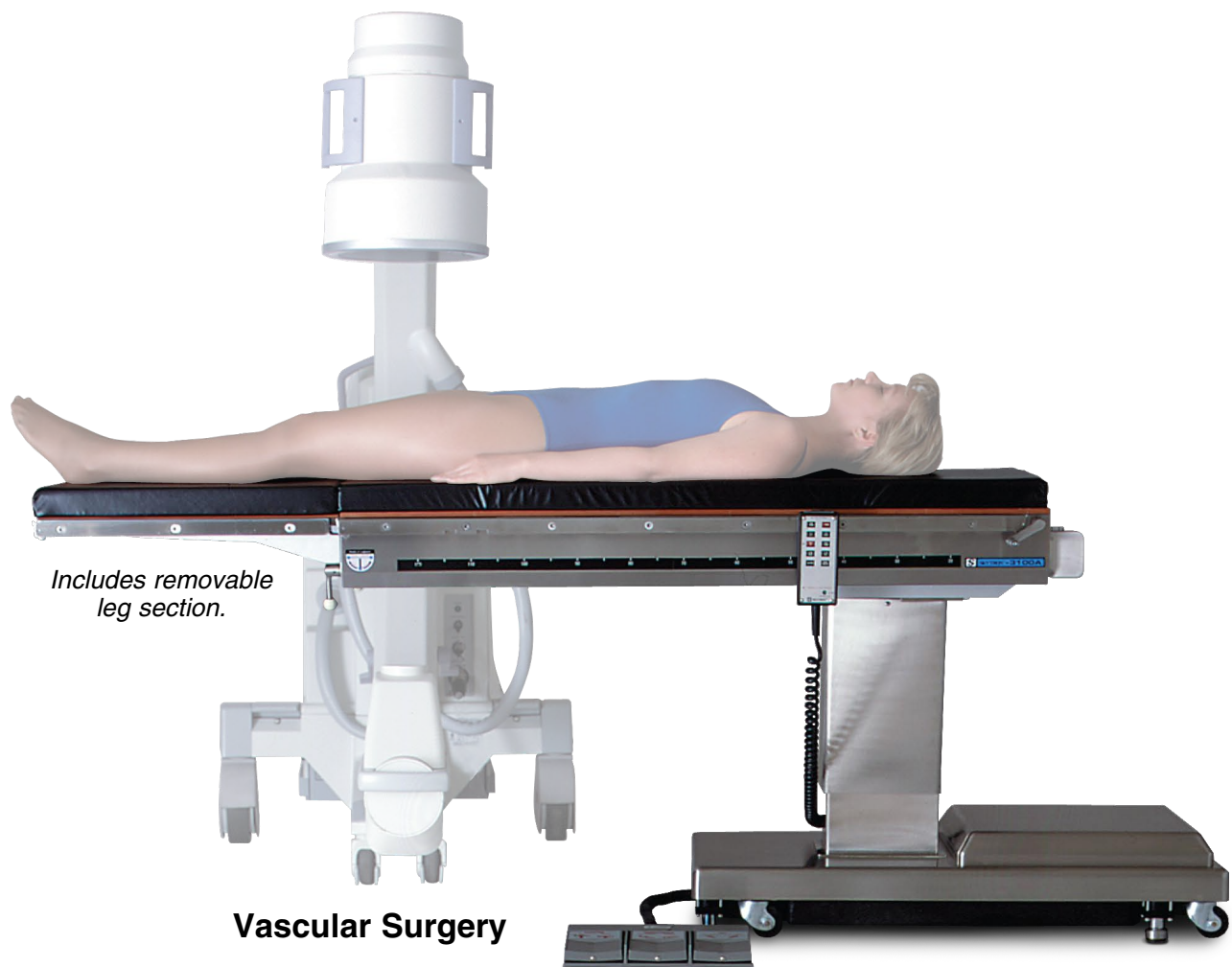

C-Arm is provided courtesy of OEC Medical Systems.

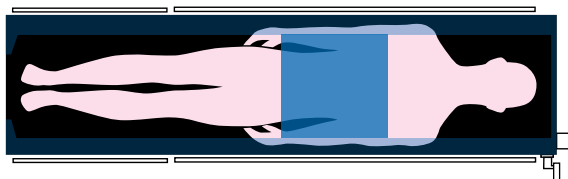


Normal Position



23-1/2 inch
Top Slide



40 inch Max.
29 inch Min.
Table Height

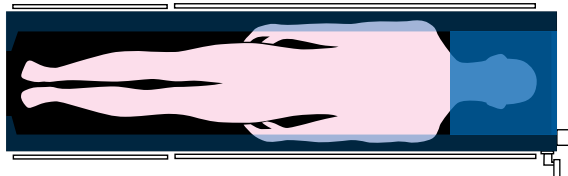

Top Slide



20° Normal and Reverse
Trendelenburg




20° Lateral Tilt
Right & Left

Positioning Capability to Compliment a Variety of Surgical Procedures



Urology
(Removable leg section)

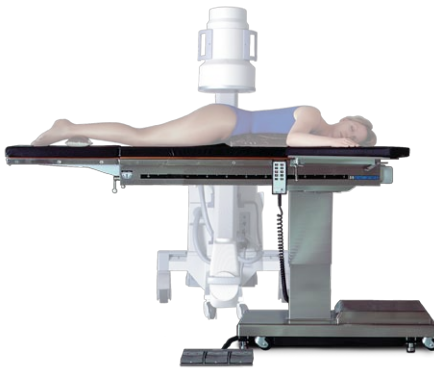

Pain Management




**500 lb. Lift and
Articulation Capacity**

One Touch Pendant Control


- Recessed buttons reduce the possibility of accidental activation
- Durable, water resistant steel housing
- Pendant control can be used from any point around the perimeter of the table
- Brakes automatically lock when any button is pushed for patient and staff safety
- Can be used in conjunction with optional foot control



**Cassette Tunnel
with Manual Crank**

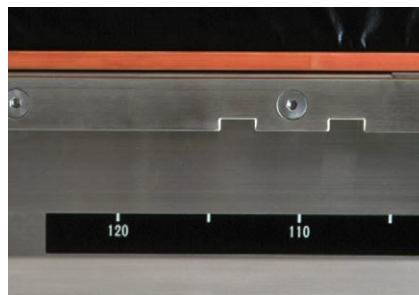

Cassette Travel Locator

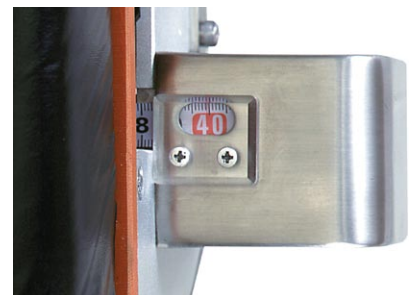

Cassette Travel Indicator

A cassette tunnel with manual crank allows convenient loading of x-ray cassettes at the head of the table, and a cassette travel indicator shows precise location of cassette.

Elite 3100A Top Slide Surgical Table Specifications

Technical Specifications:

Length	77 inches	1955mm
Width with side rails	22 inches	560mm

Adjustments via pendant control:

Height	29-40 inches	736-1016mm
Lateral tilt (left/right)	20°	
Trendelenburg	20°	
Reverse Trendelenburg	20°	
Return to level (anatomically correct return to horizontal position from articulated positions)		
Top Slide	23-1/2 inches	
Base braking (lock/unlock)		

Weight of operating table: 774 lbs. 351kg

Max. permissible load: 500 lbs. (226kg)
Lift and Articulation



Specialty top slide, electrohydraulic line powered surgical table features 500 lb. (226kg) lift & articulation capacity. Table is equipped with a 15 foot (457cm) power cord. Four (4) conductive casters. Four point braking system with push button operation, including manual override for safety and stability. A cassette tunnel with travel indicator and easy-to-load and position tray is built into the table top carriage for full image intensification capability. Hand held pendant control includes steel housing design providing extreme durability. Pendant permits control at any point around table's perimeter.



5085 Corporate Exchange Blvd. SE • Grand Rapids, MI 49512
1-800-SKYTRON (759-8766) • email: sales@skytron.us

Local Representative

MAR-D-0021 REV0
9/10