



EZ Slide 3501C



From the Leader in Top Slide Surgical Tables

On Demand Top Slide at the Touch of a Button!



The EZ Slide 3501C Top Slide surgical table is designed to meet O.R. suite demands of the 21st Century, including an impressive range of full body imaging and articulating table functions as required for C-arm imaging, without multiple accessories or having to reposition the patient. One touch of a button does it all!



7-5/8 inch Slide to Head



13-1/8 inch Slide to Foot



44-7/8 inch Max.
26 inch Min.
Table Height



21 inch
Top Slide



90° Reflex



40° Flex



98° Leg Down



90° Back Up



40° Back Down



25° Trendelenburg




20° Reverse
Trendelenburg



20° Lateral Tilt
Right & Left

EZ Slide 3501C Imaging Versatility

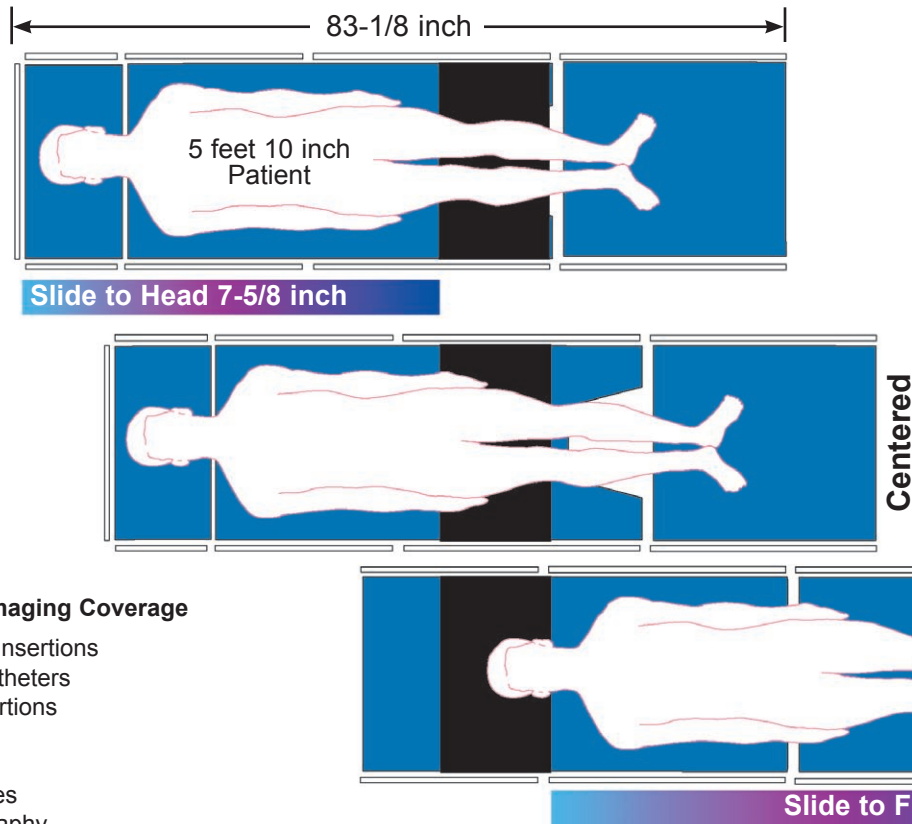

Upper Body Imaging Coverage

- Pacemaker Insertions
- Hickman Catheters
- Broviac Insertions
- Shunts
- IVC Filters
- Lung Biopsies
- Cholangiography
- Cardiovascular
- Thoracic
- Kidney
- Gall Bladder
- Laparoscopy
- Total Hip Revisions
- Ophthalmic
- Orthopedic
- Upper Limb Arthroscopy
- Neurosurgery
- Spinal
- Implantable Defibrillation
- Endoscopy

Lower Body Imaging Coverage

- Spinals
- Femoral Popliteal Bypass
- Laser Angioplasty
- Endovascular
- Orthopedics
- Lower Limb, Non-traction Orthopedics
- Laparoscopy
- Upper Limb Arthroscopy
- Cystoscopy
- Endourology
- Gynecology
- Abdominal
- Nissen Fundoplication



Durable, Back-Lit Pendant Control

- Simple One Touch Operation
- Back-Lit for High Visibility
- Audible Alarm for Kidney Bridge, Back Up Articulation, Top Slide Obstruction
- Anatomically Correct Return to Level
- Top Slide Return to Center Indicator
- Recessed Buttons to Eliminate Accidental Activation
- Exceptionally Durable (including superior steel housing design), Water Resistant Pendant Control
- Can be used in conjunction with optional Three Function Foot Control

EZ Slide 3501C Provides Optimal Imaging



**700 lbs. Lift
and 600 lbs.
Articulating Capacity**



Battery Life Indicator



Adjustable Head
Section with Quick
Release Bar for
Emergencies



Full Length X-Ray (Bucky)
Top Accessory for Image
Capture



Secondary Function
Switches for Backup to
Primary Control



Power Kidney Bridge
for Safe and Easy
Positioning

Coverage and Top Slide Range!



***Exclusive 21 inch
Top Slide***



Optional 40 inch Carbon Fiber Extension
for Specialty Imaging Requirements



Durable, Back-lit
Pendant Control



Manual Emergency
Brake Unlock



Self Leveling Brakes
for Stability on Uneven
Floors



Unique Table
Accessories for
Safe and Improved
Positioning

EZ Slide 3501C

Positioning Guidelines



Gall Bladder / Thoracic



Vascular / Endovascular



Gastro / Intestinal



***GYN / Endourology /
Cystoscopy***



Kidney / Thoracic




Ophthalmic / ENT



Neurosurgery



Spinal

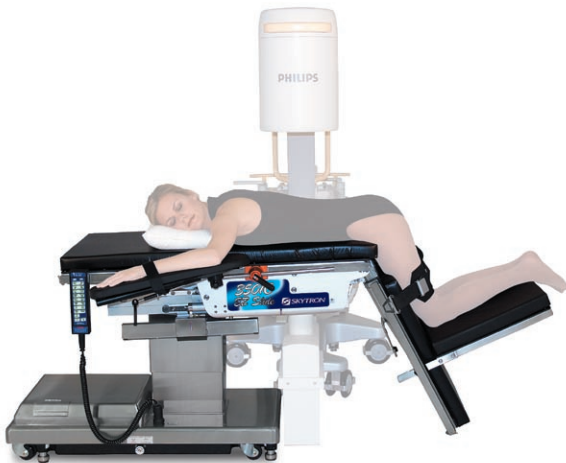

Lumbar



Shoulder / Arthroscopy



Bariatric Surgery



**Lower Extremity / Endovascular /
Podiatry**

EZ Slide 3501C

Surgical Table Specifications

Technical Specifications:

Length w/o head rest	69 inches	1750mm
Length with head rest	83-1/8 inches	2111mm
Width (back/seat plate)	19-3/4 inches	500mm
Width with side rails (back/seat plate)	21-7/8 inches	555mm



Adjustments via pendant control:

Height	26 to 44-7/8 inches	660-1140mm
Lateral tilt (left/right)	20°	
Trendelenburg	25°	
Reverse Trendelenburg	20°	
Power Kidney Bridge	6 inch Lift	150mm Lift
Back up/down	90°/40°	
Leg up/down	0°/98°	
Return to level (anatomically correct return to horizontal position from articulated positions)		
Top slide/head	7-5/8 inches	195mm
Top slide/foot	13-1/8 inches	335mm
Flex/Reflex	40°/90°	
Base braking (lock/unlock)		

Adjustments via manual control:

Head Section	60° up/90° down	
Weight of operating table:	827 lbs.	376kg
Max. permissible load:	700 lbs. (315kg) Lift	
	600 lbs. (270kg) Articulation	

Full function, electrohydraulic battery powered surgical tables feature 21 inch (530mm) longitudinal top slide, 700 lbs. (315kg) lift and 600 lbs. (270kg) articulation capacity. Table is equipped with a removable 15 foot (4.5m) power cord for battery charging. Four (4) conductive casters and four point braking system with push button operation, including manual override for safety and stability. Hand held, back-lit pendant control includes steel housing design providing extreme durability. Pendant permits control at any point around table perimeter.



5085 Corporate Exchange Blvd. SE • Grand Rapids, MI 49512
1-800-SKYTRON (759-8766) • email: sales@skytron.us

Local Representative

MAR-D-0026 REV1
1/12