



# AURORA FOUR



 **SKYTRON<sup>®</sup>**  
TRUSTED PARTNER | PROVEN SOLUTIONS

# Putting the surgical team in control of focus & illumination.

After spending decades lighting the way for surgical teams around the world, Skytron has combined all of the features important to the operating room in the Aurora Four light. This next-generation surgical light puts enhanced clarity and focus in the clinician's hand – with only one twist of the center sterile handle.

## BETTER, BRIGHTER LIGHT

Now featuring a larger spot size and intensity up to 160,000 lux, the Aurora Four surgical lights bring crisp, deep-cavity illumination to any procedure. With unique reflector technology, the lights dramatically minimize shadows that obscure visibility of the surgical site.

## EFFORTLESS OPERATION

Featuring focus control from the center sterile handle, seeing precise detail no longer means repositioning the lighthead, adjusting the spot size or leaving the sterile field. The focal point adjusts according to the surgical team's preference – with only one twist of the handle.

## MORE OPTIONS

Skytron puts the clinician in control with the option to choose between a warmer or cooler color temperature based on the procedure type. Integration with booms, flat-screen displays and HD cameras make the Aurora Four surgical light a versatile solution.

DIRECT  
FOCUS  
CONTROL

BRILLIANT  
& CRISP  
ILLUMINATION

SELECTABLE  
COLOR  
TEMPERATURE



## Direct Focus Control

### THE SURGICAL LIGHT WITH TRUE FOCUS ADJUSTMENT

When precise detail must be visualized at a greater working distance, the clinician can simply twist the center sterile handle for increased focus. The Aurora Four surgical lights put the surgical team in control while keeping the operating room hassle free.



One twist of the handle moves the focal point up or down.

### DIFFERENTIATING FEATURES

- Center sterile handle adjusts intensity, focus and positioning
- Vertically segmented reflector design for shadow reduction
- 24" or 30" diameter lightheads
- Maximum intensity of 160,000 lux
- Selectable color temperature of 4100K or 4500K
- 360° rotation on all axes
- Camera-ready options available





## Thoughtfully designed for today's surgical teams


Skytron knows that operating room services are complex. That is why the newest generation of surgical lights is designed to take unnecessary work out of each procedure.

### HASSLE-FREE FEATURES INCLUDE

- Aurora Four lights can be placed further away from the patient, opening up the line-of-sight and minimizing physical contact with the fixture
- 110° vertical travel allows downward positioning for hip, OB-GYN and urology procedures
- 360° rotation on all axes and lightweight construction keep positioning simple
- Drift-free stability keeps lights exactly where needed
- Reliable performance and lower cost of ownership bring peace of mind



Combination of fixed and variable focus LEDs produce a stunningly brilliant light

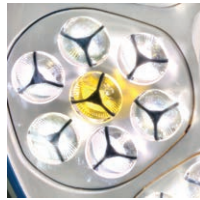

Unique vertically segmented reflectors create a crisp, homogenous spot



Wall intensity control offers five adjustable settings and color temperature selection



Color temperature of 4100K illuminates deep tissue with soft, white light



Color temperature of 4500K illuminates orthopedic procedures with bright, white light



Optional Precision HD camera system



## BRILLIANT AND CRISP ILLUMINATION

- Focus adjustment extends the working distance
- Unique upward-facing LEDs are suspended in cradles and keep the surgical field cool and comfortable
- Large spot size fills cavity with seamless light
- Maximum intensity up to 160,000 lux
- CRI of 96





## More Options for More Environments







Seamlessly integrate Aurora Four surgical lights with booms, flat-screen displays or independent camera arms. For specialty rooms or clinics, a mobile stand option is available.

### ALSO USE AURORA 4 LIGHTS IN:

- General and Hybrid Operating Rooms
- C-Section Suite
- Cath Lab
- Emergency Department
- Intensive Care Unit
- Outpatient Surgery Centers
- And more

With numerous options and configurations, we have an Aurora Four surgical light that is right for any clinical environment. From hybrid to general operating rooms, contact your local representative to find the right solution.



Model	Configuration	Lighthead Diameter	Camera-Ready Lighthead	Flat-Screen Display Arms Available	Camera Arm Available	Extended Arm Available	Laminar Flow Hub Available	Boom Integration Available	Mobile Stand Available
									
<b>AUA5</b>	Single	24"	No	1 or 2	Yes	Yes	Yes	Yes	Yes
<b>AUA55</b>	Dual	24"	No	1 or 2	Yes	Yes	Yes	Yes	No
<b>AUA555</b>	Triple	24"	No	1	Yes	No	No	No	No
<b>AUA5TV</b>	Single	24"	Yes	1 or 2	No	Yes	Yes	Yes	No
<b>AUA55TV</b>	Dual	24"	Yes	1 or 2	No	Yes	Yes	Yes	No
<b>AUA555TV</b>	Triple	24"	Yes	1	No	No	No	No	No
<b>AUA7</b>	Single	30"	No	1 or 2	Yes	Yes	Yes	Yes	No
<b>AUA75</b>	Dual	30", 24"	No	1 or 2	Yes	Yes	Yes	Yes	No
<b>AUA755</b>	Triple	24", 30", 24"	No	1	Yes	No	No	No	No
<b>AUA7TV</b>	Single	30"	Yes	1 or 2	No	Yes	Yes	Yes	No
<b>AUA7TV5</b>	Dual	30", 24"	Yes (30")	1 or 2	No	Yes	Yes	Yes	No
<b>AUA75TV</b>	Dual	30", 24"	Yes (24")	1 or 2	No	Yes	Yes	Yes	No
<b>AUA57TV5</b>	Triple	24", 30", 24"	Yes (30")	1	No	No	No	No	No
<b>AUA575TV</b>	Triple	24", 30", 24"	Yes (24")	1	No	No	No	No	No

# AURORA FOUR

INNOVATIVE SOLUTIONS FOR THE OR, ER AND EVERY ROOM IN BETWEEN



 **SKYTRON®**  
**TRUSTED PARTNER | PROVEN SOLUTIONS**



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.us](mailto:info@skytron.us)  
[www.skytron.us](http://www.skytron.us)

**About Skytron** Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at [www.skytron.us](http://www.skytron.us).