



## *Essentia 1700*



**Quality, Value,  
Durability & Performance**

# Economically Priced to Deliver Quality, Value and Durability

*Skytron's Model 1700 Surgical Table is engineered to deliver quality, value and long term durability, yet is economically priced. The 1700 general purpose surgical table is built to the same rigorous standards as our more expensive Skytron speciality surgical tables, equipped to provide full body imaging windows with 500 lb. capacity and a menu of essential table features.*




- 500 lb. Capacity
- Full Body Imaging
- 180° Top Rotation
- One-Touch Pendant Control
- Line Power
- Manual Emergency Brake Release



**180° Top Rotation provides Full Body Imaging & Uninterrupted C-Arm approach as well as clearing the base for seated Surgeons.**

**One Touch Pendant Control**

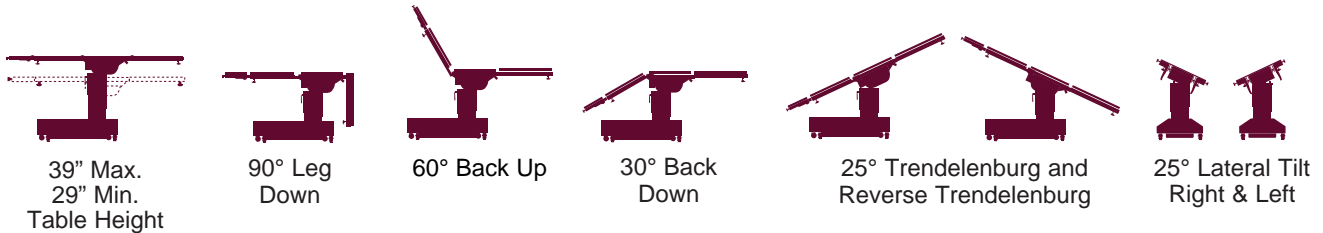

# Essentia 1700 Positioning Guidelines



Upper Body Imaging



ENT



Shoulder Arthroscopy



Knee Arthroscopy



OB/GYN & C-Section



Rectal/Hemorrhoidectomies

# Essentia 1700 Surgical Table Specifications

## Technical Specifications:

Length w/o head rest	65.5 in.	1664mm
Length with head rest	76 in.	1930mm
Width	19.75 in.	500mm
Width with side rails	22 in.	560mm



## Adjustments via hand control:

Height	29-39 in.	737-990mm
Lateral tilt (left/right)	25°	
Trendelenburg	25°	
Reverse Trendelenburg	25°	
Back up/down	60°/30°	
Leg down	90°	
Base braking (lock/unlock)		
Weight of operating table	530 lb. (240kg)	
Max. permissible load	500 lb. (227kg)	Capacity

Electrohydraulic line powered surgical table features 180° top rotation, 500 lb. (227kg) capacity. Table is equipped with a removable 15 ft. (457cm) power cord. Four (4) conductive casters, four point braking system with push button operation, including manual override, for safety and stability. Hand held pendant permits control at any point around table's perimeter.



Local Representative