

I-Series Washer-Disinfectors



***High Capacity Cleaning & Disinfection
for Surgical Instruments and Utensils***

High Capacity, Quick Processing Times & Flexible Cycle Options

The Skytron I-Series Washer-Disinfector combines high capacity cleaning with quick processing times and flexible cycle options while maximizing throughput to meet the demands of today's healthcare facilities.




I5 Single Door or Pass Thru Washer-Disinfector

I5 Washer Disinfector utilizes a highly visible, easy to use keypad control

Added Benefits:

- Volumetric Flow Meters Standard for Chemical Pumps
- Efficient Utilization of Water Channels and Spray Arms for Drying Process

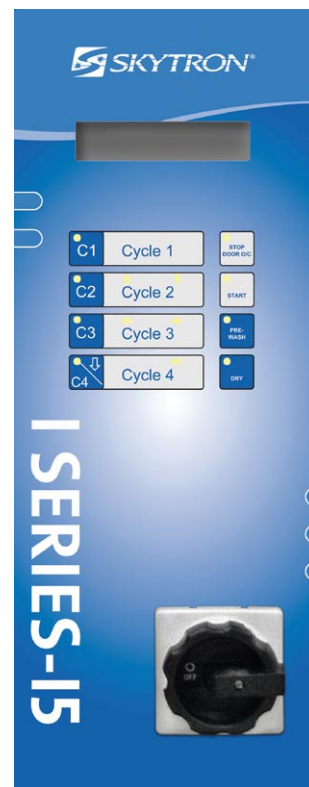

Multiple Rack Configurations to Meet Processing Needs

Easy to Use, Quiet & Efficient Operation

The Skytron I-Series Washer-Disinfectors are easy to use and engineered for quiet operation. System design features greater processing options by allowing 16 different programmed cycles. Surgical and cannulated instruments, anesthesia items, basins, pans and utensils can be processed thoroughly and efficiently.

Skytron I-Series Washer-Disinfector

Key Features Overview

I5 Key Features

- Type 304 Stainless Steel Construction
- Available in Single (S) and Pass Thru (PT) Configurations
- Power Drop Down Vertical Doors
- Full Viewing Window(s) with Low-E Glass
- Interior Halogen Light
- Integrated Dryer Package
- Two Standard Chemical Injection Pumps with Volumetric Flow Meters
- Sixteen Pre-Programmed Customizable Cycles
- Cold Water Pre-Wash
- Hot Water Detergent Wash
- Temperature Guarantees for Each Step of Every Cycle
- Thermal Disinfection Rinse (180-194° F)
- Forced Hot Air Dry Utilizing Spray Arms
- Clean Side Cycle Data Printer (optional)
- Multiple Basket Configurations to Accommodate Surgical Instrumentation, Anesthesia Equipment and Miscellaneous Equipment
- Low Water Consumption



I5 Washer Disinfector

I3 Key Features

- Type 304 Stainless Steel Construction
- Full Viewing Window with Low-E Glass
- Integrated Dryer Package
- Sixteen Pre-Programmed Customizable Cycles
- Cold Water Pre-Wash
- Hot Water Detergent Wash
- Temperature Guarantees for Each Step of Every Cycle
- Thermal Disinfection Rinse (180-194° F)
- Forced Hot Air Dry Utilizing Spray Arms
- Cycle Data Printer (optional)
- Multiple Basket Configurations to Accommodate Surgical Instrumentation
- Low Water Consumption
- Undercounter Design



I3 Washer Disinfector




I-Series Washer-Disinfector Dimensions, Options & Material Handling Accessories

Model	Interior Cabinet Size	Exterior Cabinet Size
I3	20"W x 21"D x 21.3"H 510mm x 535mm x 540mm	31.5"W x 26.8"D x 33.5"H 800mm x 680mm x 850mm
I5	25.5"W x 27.25"D x 25.25"H 648mm x 692mm x 641mm	35.5"W x 33"D x 79"H 900mm x 840mm x 2000mm

- I5 Washer-Disinfectors are available in Single Vertical Door (S) or Double Door Pass-Thru (PT)
- I3 Washer-Disinfectors available with manual drop down door

I-Series Material Handling Accessories

- Transfer Cart
- 2 Level Processing Rack
- 3 Level Processing Rack
- 4 Level Processing Rack
- 5 Level Processing Rack (I5 only)
- Spiral Anesthesia Processing Rack
- MIS Processing Rack
- Rigid Instrument Container Rack
- Pass-Thru Return Window
- Processing Rack Inserts for Miscellaneous Items

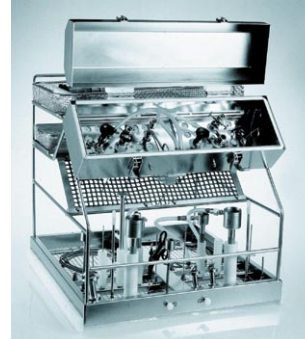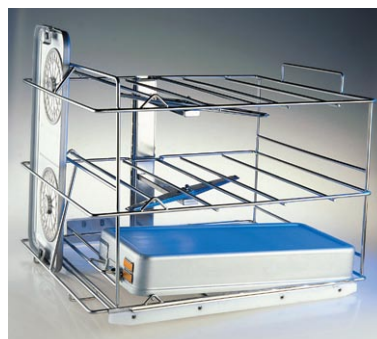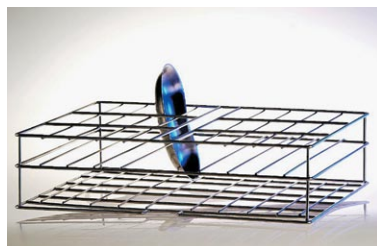

5085 Corporate Exchange Blvd. SE, Grand Rapids, MI 49512
 1-800-SKYTRON (759-8766) Fax: 616-656-2906
 e-mail: sales@skytron.us

Local Representative