



STELLAR XL



 **SKYTRON[®]**
TRUSTED PARTNER | PROVEN SOLUTIONS


High value surgical lighting with high-tech capabilities.

Well-suited for any surgical application, the Stellar XL light delivers advanced focusability and optimal shadow control to the operating room. The Stellar XL offers the ideal balance of features and value – whether updating outdated technology or outfitting a new space.

HANDLE CONTROLLED FOCUS CAPABILITY

Featuring a 27" range of total focusable depth, the Stellar XL's focal point can be easily fine-tuned with the center sterile handle – without moving the lighthead itself


- Set and forget the lighthead at a specific height
- Extend the working distance without moving the lighthead
- Reduce risk of contact with the lighthead or impeding line of sight to staff or displays

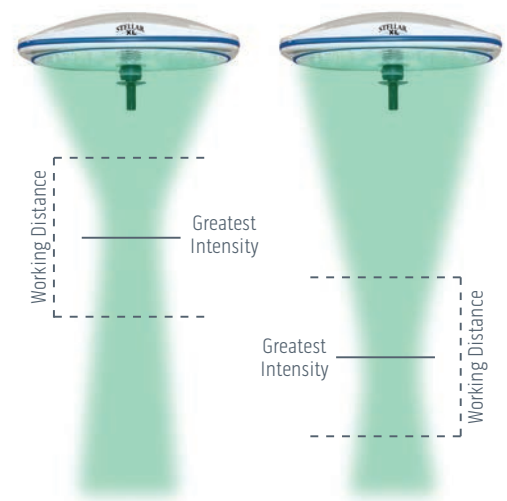

SUSPENDED LED SYSTEM

Each of the high-intensity LED is suspended upward-facing in its cradle, allowing the vertically segmented reflector (VSRD) to redirect overlapping rays down to the surgical field. The result is 2,000+ beams of light that converge into a continuous column, leading to:

- A crisp, homogeneous spot
- Excellent shadow reduction
- Hidden LEDs that eliminate glare
- Improved thermal control
- A cool, comfortable working environment

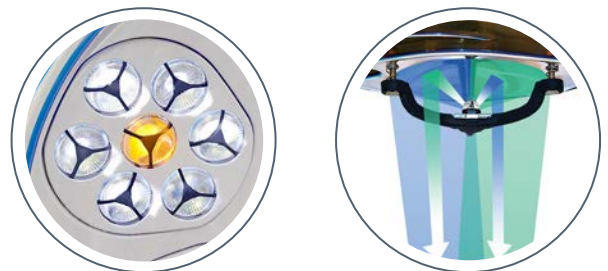

TRUE
FOCUS
ADJUSTMENT

VERTICALLY
SEGMENTED
REFLECTOR
DESIGN

SUSPENDED
LED
SYSTEM



Ideal For Any Application

Available in single, dual or triple lighthead configurations and fully compatible with Skytron booms, the Stellar XL is adaptable to many clinical settings.

PERFECT FOR

- Upgrading outdated technology
- Outpatient Surgery Centers
- General Operating Rooms
- C-Section Suites
- Labor and Delivery Suites
- Emergency Departments
- Intensive Care Units
- Procedure Rooms



Also available on a mobile stand

DIFFERENTIATING FEATURES

- Maximum intensity up to 160,000 lux
- Natural warm color temperature provides white light for high intensity applications
- 27" range of total focusable depth
- 24" diameter lighthead
- 360° rotation on all axes
- Lightweight polymer housing permits easy maneuverability

STELLAR XL

INNOVATIVE SOLUTIONS FOR THE OR, ER AND EVERY ROOM IN BETWEEN



The literature has been developed to provide healthcare professionals with general product information. Please note: not all products are available in all countries. Product labeling and associated claims may differ based on regulatory requirements. For product availability or any further information, please contact your local Skytron representative.

 **SKYTRON®**
TRUSTED PARTNER | PROVEN SOLUTIONS

DISTRIBUTED BY
MEDITEK

Here for L.I.F.E.

1.800.567.8400 www.meditek.ca
1810 Ellice Avenue Winnipeg, MB R3H 0B7

About Skytron Skytron is a privately held company specializing in capital equipment that drives efficiency in today's healthcare facilities. Our surgical tables, surgical & exam lighting, booms, video integration systems and other products help hospitals around the world work better. To learn more, visit our website at www.skytron.us.